


ANZAI

since 1976

Anzai Breath Learning Equipment

ABLE

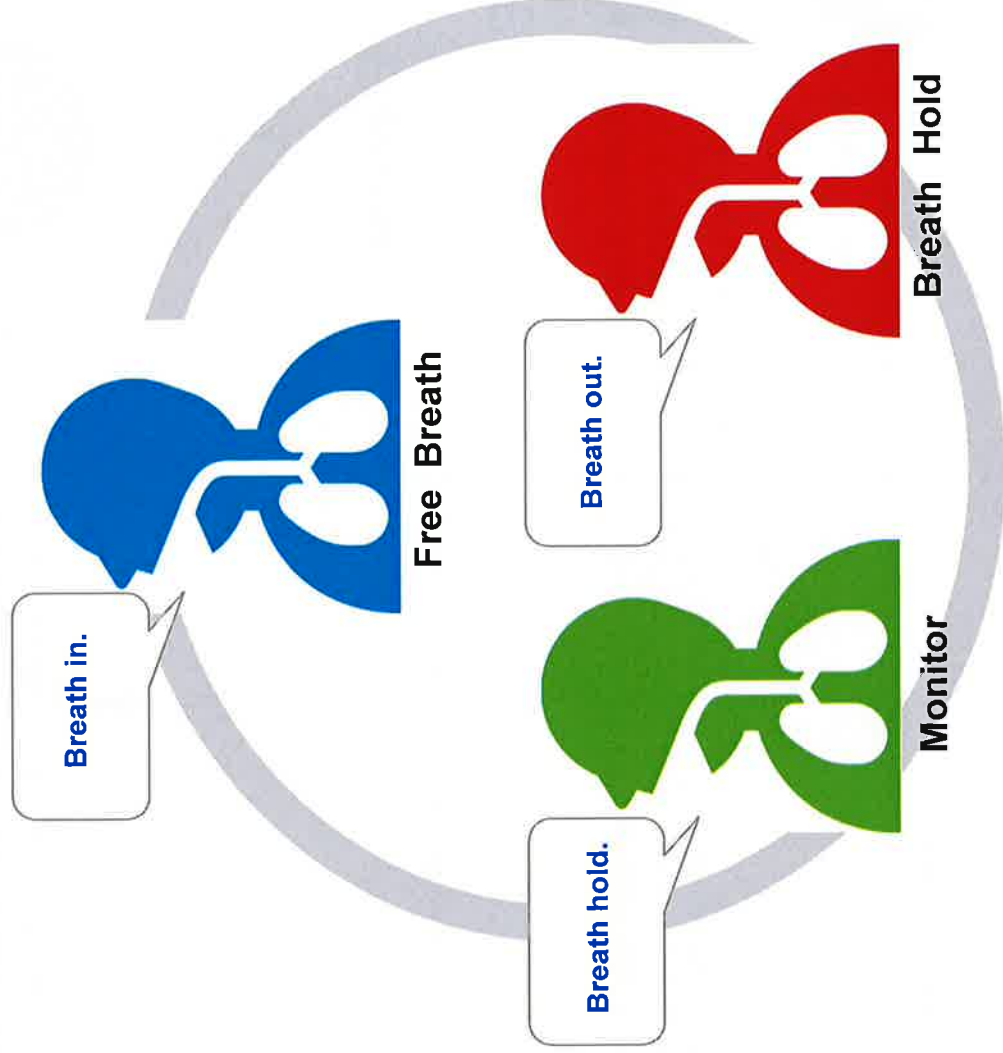


ABLE is a breath stabilization monitor intended for breath coaching (free breath, breath hold) by utilizing a Tablet and a Goggle.

Visual (image) information and acoustic (voice) information are adopted for breath coaching methods.

It enables breath coaching without the limitation of time and place such as clinical sites where patient respiration is relied on and breathing exercise at home.

ABLE is capable of displaying the respiratory waveform from the Anzai Respiratory Gating System AZ-733VI. This combination use with ABLE will be expected to realize more effective diagnosis and treatment.

Breath Coaching Screen

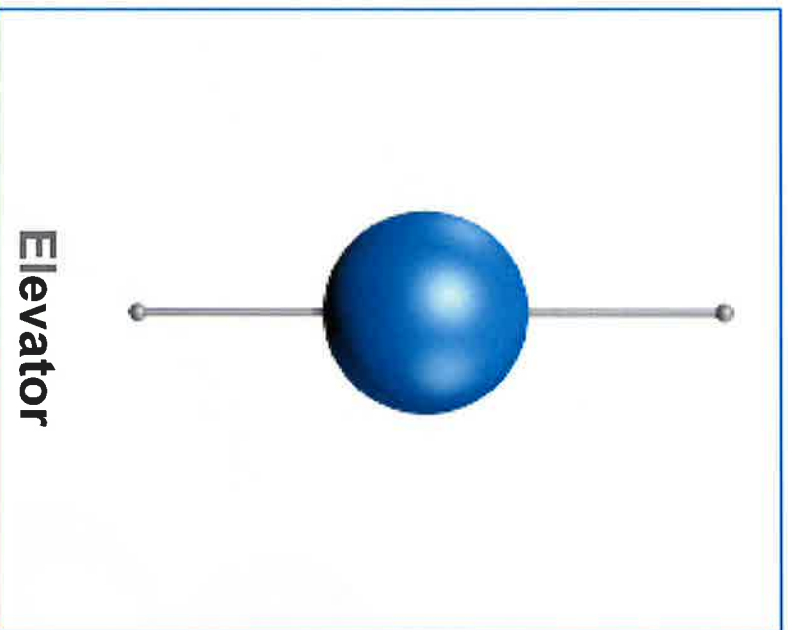

Elevator

Breath coaching video is provided to the Tablet. * 1
For example, “Elevator” coaches breath with ball movement. When ball is moving upward, it tells the timing of inhalation. When moving downward, then exhalation.

By providing a simple video, anyone from children to elderly persons can perform breath exercise.

Eight languages are supported for both character and voice display: Japanese, English, Russian, Spanish, French, German, Portuguese and Chinese). * 1

* 1 Breath coaching video and supporting language is to be increased accordingly

Component




■ Tablet *2 2 unit
ABLE application is installed in the Nexus6 (platform: Android OS, screen size: 6.0 inches)



■ Goggle 1 unit
This is applied to patients.



■ Wireless transceiver 1 unit
This is used to display the respiratory waveform obtained by the Respiratory Gating System AZ-733VI on the ABLE. It supports Bluetooth.

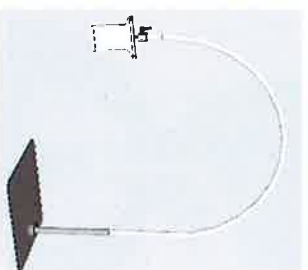

■ Tablet Fixture 1 unit
This is a dedicated flexible arm stand for the Tablet.

*2 Application for iPhone is under preparation.

This product subjects to change without notice.



Anzai Medical Co., Ltd.

Add : 3-6-25 Nishi-Shinagawa, Shinagawa-ku Tokyo, 141-0033, Japan

Tel : +81-3-3779-1611 Fax: +81-3-3779-6606

Email : info@anzai-med.co.jp Web : www.anzai-med.co.jp